


# ASMR: CADD Medicine Bag Change



## STEP 1 – SET UP

Clean the area where your supplies will be placed using alcohol or a disinfectant wipe and gather your supplies.

- (1) Alcohol-chlorhexidine Swabs
- (1) New Medicine Bag
- (4) AA Batteries if due to be changed

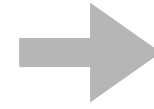


## STEP 2 – CLEAN HANDS

Wash your hands with soap for 30 seconds. Dry with paper towel. Turn tap off with paper towel.

Got Questions?  
**Business Hours:**  
Call or text: 604-862-6824  
Email: [pivotclinic@cw.bc.ca](mailto:pivotclinic@cw.bc.ca)  
**After Hours:**  
Call: 604-875-2445 and say "IV Team" when prompted (weekdays after 5pm and weekends)

# A S

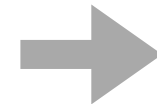

# M




# R

## STEP 3 – Alcohol Swab

1. **Open** alcohol-chlorhexidine swab
2. **Check** New Medicine Bag\*
3. **Remove** pump and old medicine bag from pouch/velcros.
4. **Pause** CADD pump:
  - a. press start/stop button
  - b. press yes to pause
5. **Sanitize** hands
6. **Scrub** CADD tubing spike with alcohol-chlorhexidine swab where it is inserted into the old medicine bag for 30 seconds\*\*
7. **Allow to dry** for 1 minute

swab

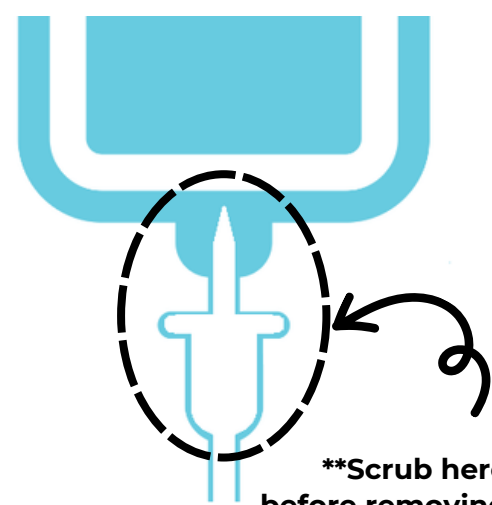

\*\*Scrub here before removing spike from old medicine bag

### \*\*SCRUB AREA:

Scrub between the used medication bag and the tubing. Scrub the connecting area and the inch above and below the connecting area for 30 seconds. Allow to dry for 1 minute.

## STEP 4 – Medicine

1. **Remove tab** from new medicine bag
2. **Remove CADD tubing spike** from old medicine bag - be careful not to touch the spike
3. **Insert spike** into new medicine bag



### \*Checking New Medicine Bag:

- Name
- Antibiotic
- Dose
- Expiry
- No bubbles larger than a dime

## STEP 5 – Reset Reservoir

1. **Reset the reservoir volume on the CADD pump:**
  - a. press tasks
  - b. press reset reservoir volume
  - c. press yes
2. **Restart CADD pump:**
  - a. press start/stop button
  - b. press yes to start the infusion
3. **CADD Pouch:**
  - a. Secure CADD pump and medicine bag into carry pouch
  - b. Coil extra tubing and secure in pouch using built-in velcro straps



## STEP 6 - Battery Change (if due)

Changing Batteries (every 2-3 days):

1. **Pause & turn off** CADD pump
2. **Open** battery compartment & **replace** batteries
3. **Close** battery compartment & **turn on** CADD pump
4. **Restart** CADD pump:
  - a. press start/stop button
  - b. press yes to start the infusion

